


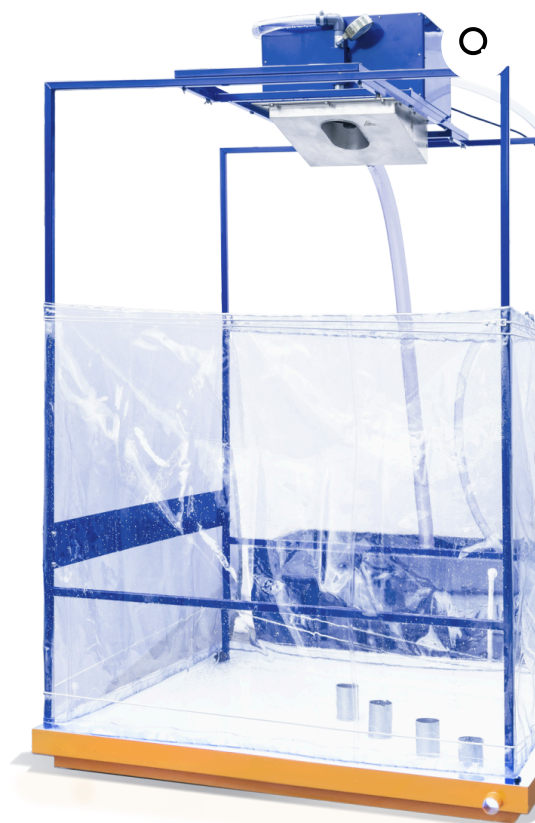

FEL3 – Rainfall Simulator

The FEL3 Rainfall Simulator can be used in the laboratory or in the field for a wide range of research from studies of infiltration under sprinkler irrigation to estimating soil loss in high intensity tropical storms.

PRODUCT CODE: FEL3

Categories: Agricultural, Educational, Environmental, FE Series - Irrigation Water Management, FEL3 - Rainfall Simulator, Geology & Geography, S Series - Hydraulics & Hydrology

Tags: APERTURE ADJUSTMENT , Control Valve , Erodibility Of Soils , FE Series , FEL3 , FIELD AND LAB TEST PLOTS , Flow Meter , High Intensity Tropical Storms , In The Field , Irrigation Water Management , Kinetic Energy , Laboratory , RAIN GAUGES , Rainfall Simulator , Simulated Rainfall , Soil Loss , Splash Erosion , Spray Head , Spray Nozzle , Sprinkler Irrigation , TILTING STAND

DESCRIPTION

The FEL3 Rainfall Simulator can be used in the laboratory or in the field for a wide range of research from studies of infiltration under sprinkler irrigation to estimating soil loss in high intensity tropical storms.

Erodibility of soils can be studied in the laboratory and the influence of crop cover on the effect of rainfall can also be investigated.

It is ideal for investigating the relationship between rainfall and soil erosion, the nature of soil erosion potential on different soil types and identifying methods by which erosion may be prevented.

The Armfield Rainfall Simulator consists of a metal frame that supports the spray head assembly. The unit may be placed directly on the ground for field studies or used with its accessory tray for laboratory experiments. In use, water is pumped from the holding tank via a control valve to the spray nozzle.

Water from the spray nozzle is interrupted by a horizontal rotating disc driven by a variable speed motor. Adjustable width apertures in the disc allow some water to pass and the remainder is recirculated via the tank. For use in field or laboratory electrical and water supplies are required.

TECHNICAL SPECIFICATIONS

FEATURES & BENEFITS

Demonstration Capabilities

- Investigate the relationship between rainfall intensity and soil erosion
- Investigate the nature of soil erosion potential on different soil types
- Investigate the methods by which soil erosion may be prevented
- Determine the kinetic energy of simulated rainfall at various rainfall intensities
- Obtain a relationship between splash erosion and intensity of simulated rainfall

Ordering specification

Rainfall Simulator comprising:

- Spray head assembly supported by metal stand
- Storage tank from which water is pumped to spray head via flow meter and control valve
- Small square test plots
- Field test plots
- Tilting stand for test plots
- Sample vessels
- Rain gauges
- Complete with instruction manual

Requirements

Electrical supply:

FEL3-A: 220-240V/1Ph/50Hz

FEL3-B: 120V/1Ph/60Hz

FEL3-G: 220-240V/1Ph/60Hz

Cold water supply

Shipping Specifications

PACKED AND CRATED SHIPPING SPECIFICATIONS

Volume: 3.9m³

Gross Weight: 576kg

Dimensions

Length: 2.65m

Width: 1.60m


Height: 2.70m

Order Codes

- **FEL3-A:** 220-240V/1Ph/50Hz
- **FEL3-B:** 120V/1Ph/60Hz
- **FEL3-G:** 220-240V/1Ph/60Hz



YOU MAY ALSO LIKE...

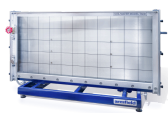

**S12-MKII-50 –
ADVANCED
ENVIRONMENT
AL HYDROLOGY
SYSTEM**

Read more

**FEL2 –
SOIL/WATER
MODEL TANK**

Read more

**S10 – RAINFALL
HYDROGRAPHS**

Read more

Armfield Limited. Unit 10 Headlands Business Park, Salisbury Road, Ringwood, Hampshire BH24 3PB.
England. Tel: +44(0) 1425 478781

